

KIZO LAB



BioGro


LED BRUSH

multi-function
hair serum applicator

USER GUIDE

**KIZO
LAB**

PARTS + CONTROL PANEL

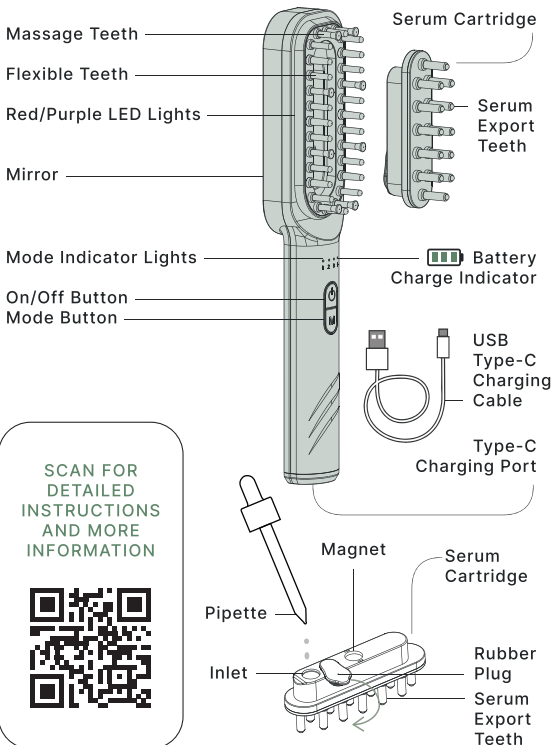

POWER UP

QUICK START

2-step treatment with **Kizo Lab Follicle Activator Hair Serum**, LED red light, LED purple light, and high intensity massage. **Before using BioGro LED Brush, boost results by applying hair serum directly onto scalp.**

STEP ONE

Fill cartridge with **Follicle Activator Hair Serum**. Press and hold power button for 3-4 seconds to turn unit on. Quickly press Mode button twice.


STEP TWO

Starting at hairline, press brush against scalp so teeth are in contact with scalp. Slowly glide brush across scalp from front to back. Repeat until all targeted areas have been covered. For best results, use twice daily.

TO FILL SERUM CARTRIDGE:


STEP ONE

Remove center Serum Cartridge from BioGro LED Brush by pulling on Export Teeth.

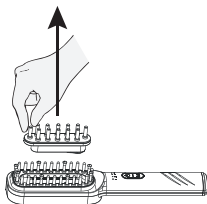

STEP TWO

Fill cartridge with **Kizo Lab Follicle Activator Hair Serum**.

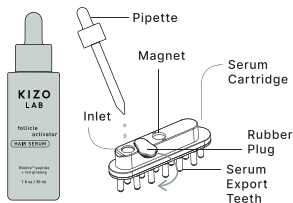

HAIR SERUM APPLICATOR + LED + MASSAGE BRUSH

- + Evenly distributes hair serum onto scalp
- + 2-speed massage brush
- + Red + purple LED
- + Automatic shutoff after 10 minutes

MADE FOR

- + people with thinning hair
- + dormant hair follicles
- + hair thickness + fullness

MODE SELECTION

MODE ONE

LED red light + serum application

MODE TWO

LED purple light + serum application +
massage vibration (LOW)

MODE THREE

LED red + purple light + serum application +
massage vibration (HIGH)

IMPORTANT SAFEGUARDS

- + Keep out of reach of children.
- + Operate this unit only according to its intended use as described in this User Guide.
- + This unit should not be used by pregnant women or individuals with heart disease, abnormal blood pressure, or cerebrovascular disease.
- + If you have questions or concerns, consult your doctor before using this unit.
- + Do not disassemble this unit.
- + Do not immerse unit in water. If unit is dropped, appears damaged, or has fallen into water, discontinue use.

SPECIFICATIONS

- + Operating current: 500-1,000 mA
- + Rated power: approx. 4W
- + Input voltage: 5V
- + Charging time: 3-4 hours
- + Product size: 222 x 45 x 33 mm
- + Battery capacity: 1,000 mAh/3.7V
- + Net weight: approx. 150 g
- + Red light wavelength: 630-660 nm
- + Purple light wavelength: 240-280 nm
- + Implementation standard: GB 4706.1-2005
GB 4706.15-2008

WARRANTY

BioGro

LED BRUSH

multi-function
hair serum applicator

LIMITED ONE YEAR WARRANTY

Kizo Lab BioGro LED Brush is warranted with normal use in accordance with User Guide instructions. Kizo Lab will replace your BioGro LED Brush free of charge for 12 months from the date of purchase if the unit is defective in materials or function. Proof of purchase is required.

EXCLUSIONS

01. This warranty does not cover any unit that has been damaged due to inappropriate usage, normal wear and tear, accident, tampering, or use with a power supply adapter other than the one provided.
02. This warranty does not cover any replacement unit if the original unit was replaced within a year from date of purchase.



CONTACT US AT 888.327.8188